

Theradiag  
More than 20 years  
experience in test  
design  
in autoimmunity

Fidis™  
Innovative multiplex  
technology with one  
of the biggest  
test menu

Immuno-trol  
Multi-parametric  
and ready-for-use  
quality controls

## One of the largest menu of multiplexed tests in autoimmunity

REFERENCE	DESIGNATION	MARKERS
MX 002	FIDIST™ THYRO*	Simultaneous detection of autoantibodies against: - Thyroperoxidase (TPO) - Thyroglobulin (TG)
MX 005	FIDIST™ dsDNA*	Quantitative detection of autoantibodies directed against double stranded DNA (dsDNA).  Simultaneous detection of autoantibodies against: - dsDNA - SS-A 60 kDa - SS-A 52 kDa (TRIM21 or Ro 52) - SS-B - Sm - Sm/RNP - Scl-70 - Jo-1 - Ribosomal Proteins - Centromere (CENP-B) - PmScl - PCNA - Histone - U1RNP
MX 117	FIDIST™ CONNECTIVE PROFILE <sup>(2)</sup>	
MX 007	FIDIST™ VASCULITIS <sup>(1)</sup>	Simultaneous detection of autoantibodies against: - Myeloperoxidase (MPO) - Proteinase 3 (PR3) - Glomerular Basement Membrane (GBM)
MX 012	FIDIST™ APS IgG	Simultaneous detection of IgG or IgM autoantibodies against: - Cardiolipin - Beta2-Glycoprotein I (β2GPI) - Prothrombin or Factor II (FII)
MX 013	FIDIST™ APS IgM	
MX 030	FIDIST™ CYTO-PANEL	Simultaneous detection of autoantibodies against: - ATPase H+/K+ of the gastric parietal cells (PCA) - Intrinsic Factor (IF) - Ribosomal Proteins - Mitochondrial M2 Proteins - LKM1 - LC1
MX 103	FIDIST™ RHEUMA-RF	Simultaneous detection of IgM class rheumatoid factor directed against Fc determinants from human or animal IgG Immunoglobulins
MX 104	FIDIST™ anti-CCP	Simultaneous detection of autoantibodies directed against 4 synthetic citrullinated peptides
MX 109	FIDIST™ CELIAC DPG IgG	Simultaneous detection of IgG or IgA autoantibodies against: - Gliadin - Deamidated Peptides of Gliadin (DPG) - tissue Transglutaminase (tTG)
MX 110	FIDIST™ CELIAC DPG IgA	

(1) FDA cleared (2) FDA cleared for CONNECTIVE 10  
All FIDIST™ kits are CE

14 rue Ambroise Croizat  
CS 90136 Croissy-Beaubourg  
F-77435 MARNE LA VALLÉE CEDEX 2

Phone +33 (0)1 64 62 10 12 – Fax +33 (0)1 64 62 09 66  
info@theradiag.com – www.theradiag.com



Watch the video !

**Theradiag**  
INNOVATION FOR BIOTHERAPIES

# Theralis

## Versatility in autoimmunity

The only system that integrates multiplex Fidis™, immunofluorescence and ELISA to answer the needs of flexibility of your core lab.

Productivity

Flexibility

Security

Pictures may differ from actual products

Read carefully the instruction for use of the product insert

Theralis - V.01/2018 - UK

**Theradiag**  
INNOVATION FOR BIOTHERAPIES

Bringing versatility to autoimmunity

Allowing the performance of up to 3 techniques simultaneously, Theralis brings you **improved productivity, flexibility** and **security** from the sample tube to the patient results.




Productivity

Full automation from primary sample to results.

- High throughput pipetting system:
- 4 ceramic coated needles,
  - 1 dedicated needle for IFA washing.

- High loading capabilities:
- up to 120 samples and 60 reagents on board,
  - up to 20 slides with 10 protocols in IFA,
  - up to 3 Fidis™ plates,
  - up to 3 ELISA microplates.

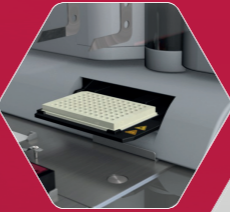


Flexibility

Theralis improves the lab workflow with:

- combination of different technologies in the same run,
- integration of ELISA microplates management,
- Fidis™ Multi-batch : association of different multiplex profiles in the same plate.

Only one instrument to process all your autoimmunity assays.

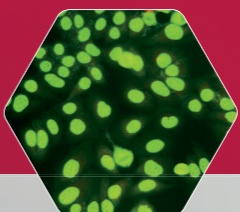

Security

- Secured loading:
- positive identification for samples and reagents,
  - lot data integration with 2D barcodes.

- Secured results:
- user login with assignable rights,
  - blocked results on invalid QC,
  - integrated Levey Jennings graphs,
  - LIS connection.

- Optimal traceability:
- lots,
  - users,
  - patient & QC results.
  - Multicriteria search engine.




The latest generation of software in autoimmunity

- One software for the whole instrument range of Theradiag in autoimmunity.
- Multitask system.
- Real time follow up of analysis.
- Total traceability: reagents, users, patients and QC.
- Integrated maintenance assistants.
- Tool for quality audit.
- Bidirectional connection to LIS.
- Drives Luminex reader and TECAN ELISA readers.

An efficient and complete solution for medium sized laboratories

Designed for lower throughput laboratories or as a backup solution, Theralis  $\mu$  gets profit of Solinium to strengthen your activity in autoimmunity.

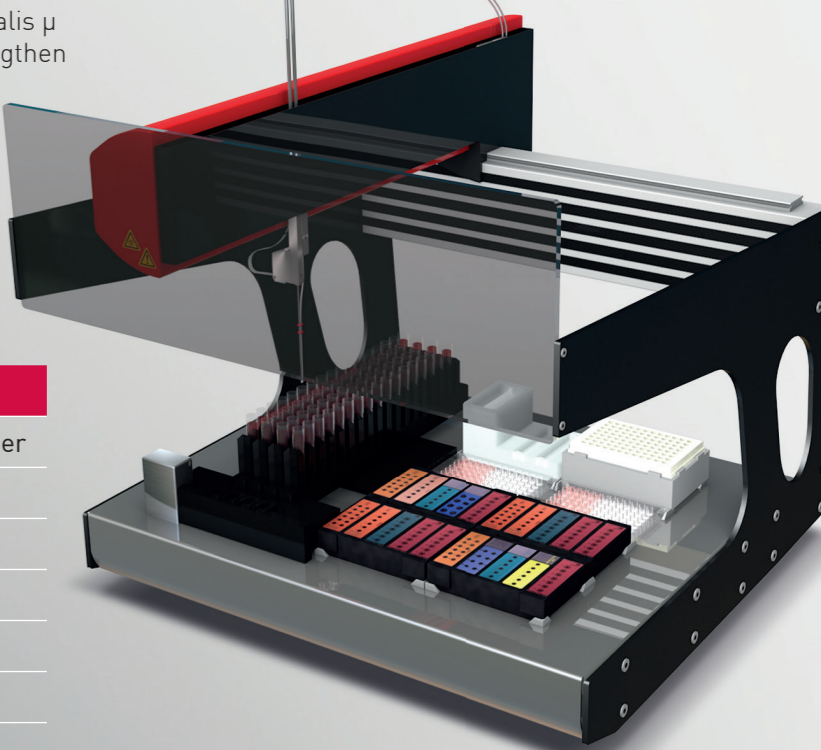

	Theralis	Theralis $\mu$
Needles	4 + 1 tip washer	1 + 1 tip washer
Samples	120	96
Reagents	60	30
Microplates	4	1
Samples positive ID	✓	✓
Reagent positive ID	✓	✗